



Dual Screen-Printed Carbon Electrode

Ref. X1110

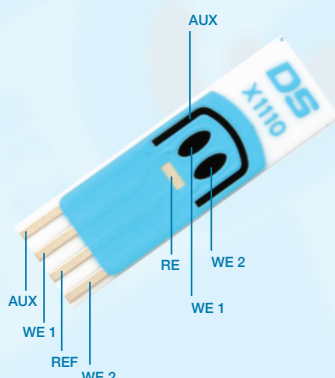

Screen-Printed Electrodes with two elliptic carbon working electrodes, one carbon counter electrode and a silver reference electrode. These electrodes may be used with DropSens portable Bipotentiostats or with any kind of bipotentiostat. They are aimed at detecting two signals simultaneously, allowing (differential) measurement of two analytes in the solution.

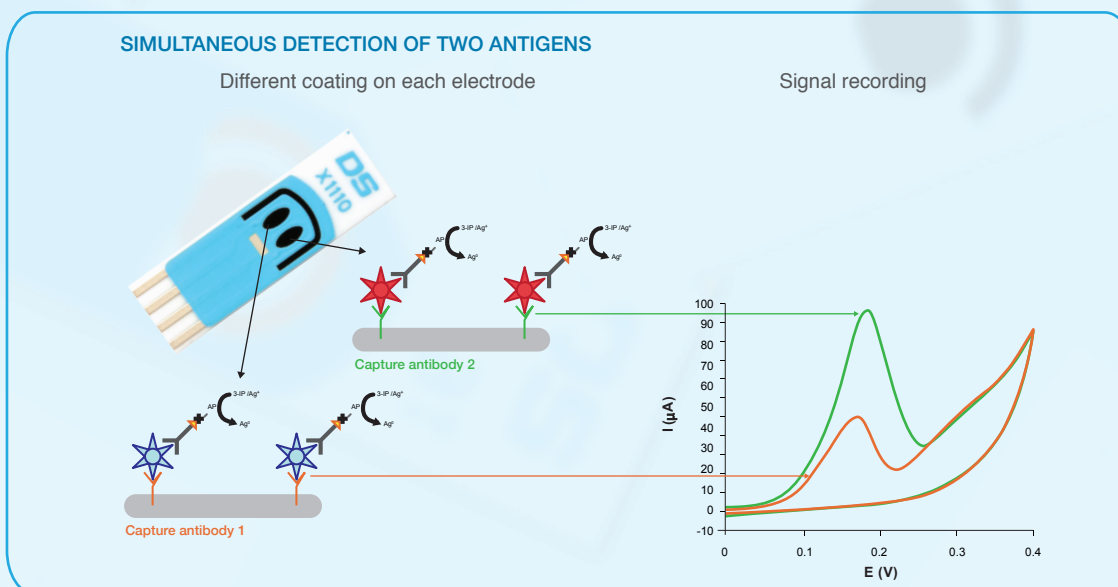

Dual Screen-Printed Carbon Electrodes are commercialised in 75 units packs. They should be stored at room temperature in a dry place.

Also specific connectors that act as an interface between the screen-printed electrode and any bipotentiostat (refs. BIDSC, BICAC) are available at [DropSens](https://www.dropsens.com).

Related products



BIDSC



BICAC



BICAST



STAT200



STAT300



STAT400

Full Catalogue



Parque Tecnológico de Asturias - Edif. CEEI. 33428 LLanera (Asturias). Spain
(+34) 985 27 76 85 - info@dropsens.com - www.dropsens.com

Contact Form

